GuideLiner™

V3 catheter

en	English/Instructions for Use.....	1
cs	Česky/Návod k použití.....	2
da	Dansk/Brugsanvisning.....	3
nl	Nederlands/Gebruiksaanwijzing.....	3
et	Eesti/Kasutusjuhised.....	4
fi	Suomi/Käyttöohjeet.....	5
fr	Français/Mode d'emploi.....	6
de	Deutsch/Gebrauchsanweisung.....	7
el	Ελληνικά/Οδηγίες χρήσης.....	8
hu	Magyar/Használati utasítás.....	9
it	Italiano/Istruzioni per l'uso.....	10
lv	Latviski/Lietošanas instrukcija.....	11
lt	Lietuvių/Naudojimo instrukcijos.....	12
no	Norsk/Bruksanvisning.....	12
pl	Polski/Instrukcja obsługi.....	13
pt	Português/Instruções de utilização.....	14
ru	Русский/Инструкция по применению.....	15
es	Español/Instrucciones de uso.....	16
sv	Svenska/Bruksanvisning.....	17
tr	Türkçe/Kullanım Talimatları.....	18
sk	Slovenčina/Návod na použitie.....	19
vi	Tiếng Việt/Hướng dẫn sử dụng.....	20
hr	Hrvatski/Upute za uporabu.....	21
sr	Srpski/Uputstvo za upotrebu.....	21
ro	Română/Instrucțiuni de utilizare.....	22
uk	Українська/Інструкції для застосування.....	23
th	ภาษาไทย/คำแนะนำสำหรับการใช้งาน.....	24

Teleflex™

Teleflex Medical LLC
3015 Carrington Mill Blvd.
Morrisville, NC 27560 USA
USA: (866) 246 6990
International: +1 919 544 8000

GuideLiner™ V3 catheter

Instructions For Use

USA CAUTION

Federal law restricts this device to sale by or on the order of a physician.

DEVICE DESCRIPTION

The GuideLiner catheter is a single lumen rapid exchange catheter offered in sizes compatible with 5F, 6F, 7F, and 8F guide catheters and may be placed over a standard length guidewire. The larger sizes of GuideLiner catheters are intended to be used within the proximal portions of the coronary vasculature to provide support and/or facilitate use of multiple interventional devices. The 150cm device has a 125cm stainless steel shaft section covered on the distal 17cm with a semi-circular polymer. The steel shaft is followed distally by a 25cm lumen section wiped with silicone.

The GuideLiner catheter has two platinum-iridium marker bands, which enable visibility while using standard fluoroscopic methods. The distal marker band is located on the distal tip. The proximal marker band is located near the collar. The device has two positioning marks located at 95cm (single mark) and 105cm (double mark) from the distal tip.


The GuideLiner catheter is delivered through a guide catheter resulting in an inner diameter that is approximately 1 French smaller than the guide catheter. The GuideLiner catheter has a proximal tab which indicates guide catheter compatibility and the resulting GuideLiner catheter inner diameter.

The GuideLiner catheter has been sterilized with ethylene oxide.

STERILE EO

SPECIFICATIONS

Model	Compatible Guide Catheter	GuideLiner Minimum I.D.	GuideLiner Tip O.D. (X)
5569 5F	≥ 5F (≥ 0.056" / 1.42mm I.D.)	0.046" (1.17mm)	0.053" (1.35mm)
5570 5.5F	≥ 6F (≥ 0.066" / 1.68mm I.D.)	0.051" (1.30mm)	0.063" (1.60mm)
5571 6F	≥ 6F (≥ 0.070" / 1.78mm I.D.)	0.056" (1.42mm)	0.067" (1.70mm)
5572 7F	≥ 7F (≥ 0.078" / 1.98mm I.D.)	0.062" (1.57mm)	0.075" (1.90mm)
5573 8F	≥ 8F (≥ 0.088" / 2.24mm I.D.)	0.071" (1.80mm)	0.085" (2.16mm)

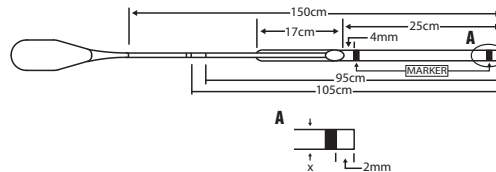INDICATIONS

GuideLiner catheters are intended to be used in conjunction with guide catheters to access discrete regions of the coronary and/or peripheral vasculature, and to facilitate placement of interventional devices.

CONTRAINDICATIONS

The GuideLiner catheters are contraindicated in vessels less than 2.5mm in diameter and in vessels in the neurovasculature or the venous system.

WARNINGS

The GuideLiner catheter is provided sterile for single use only. Reuse of single-use devices creates a potential risk of patient or user infections and may compromise device functionality, which may lead to illness or serious patient injury.

Never advance the GuideLiner catheter into a vessel without a leading guidewire as vessel damage may result.

Never advance the GuideLiner catheter into a vessel with an effective diameter less than 2.5mm. Vessel injury and/or occlusion may result. If pressure in a vessel dampens after inserting the GuideLiner catheter, withdraw the catheter until the pressure returns to normal.

Due to the size and non-tapered tip of the GuideLiner catheter, extreme care must be taken to avoid vessel occlusion and damage to the wall of the vessels through which this catheter passes.

Never advance or withdraw an intravascular device against resistance until the cause of the resistance is determined by fluoroscopy. Movement of the catheter or guidewire against resistance may result in separation of the catheter or guidewire tip, damage to the catheter, or vessel damage.

PRECAUTIONS

Do not use the GuideLiner catheter if the packaging has been damaged. A damaged package could result in a breach of sterility or device damage.

Inspect the GuideLiner catheter prior to use for any bends or kinks. Do not use a damaged catheter. Vessel damage and/or inability to advance or withdraw the catheter may occur.

The catheter lumen should be flushed with sterile, heparinized saline prior to use to prevent clot formation and to ensure the catheter is free from debris which could be introduced into the body resulting in an embolism and/or occlusion.

Precautions to prevent or reduce clotting should be taken when any catheter is used in the vascular system.

Exercise care in handling of the catheter during a procedure to reduce the possibility of accidental breakage, bending or kinking. Do not apply torque to the catheter during delivery, as catheter damage may result.

When the catheter is in the body, it should be manipulated only under fluoroscopy. Do not attempt to move the catheter without observing the resultant tip response, as catheter damage or vessel injury may occur.

Never advance the GuideLiner catheter more than 10cm beyond the tip of the guide catheter as the GuideLiner catheter may become lodged in the guide catheter making it difficult to remove.

Do not withdraw an undeployed stent back into the GuideLiner catheter when the catheter is in the body, as it may result in dislodging the stent. Instead, simultaneously pull both the GuideLiner catheter and undeployed stent back into the guide and remove together.

ADVERSE EFFECTS

Potential adverse effects that may be associated with the GuideLiner catheter include, but are not limited to, the following:

- Cardiac arrest
- Embolism
- Infection
- Myocardial infarction
- Slow-flow/Occlusion
- Stent dislodgement
- Thrombosis
- Vessel dissection
- Vessel perforation

CLINICAL PROCEDURE

The GuideLiner catheter should be used by physicians trained on the procedures for which it is intended. The techniques and procedures described do not represent ALL medically acceptable protocols, nor are they intended as a substitute for the physician's experience and judgment in treating any specific patient. All available data, including the patient's signs and symptoms and other diagnostic test results, should be considered before determining a specific treatment plan.

PACKAGE CONTAINS:

1x Catheter

OTHER ITEMS REQUIRED BUT NOT PROVIDED:

- Guide catheter with an inner diameter large enough to accommodate the specific model of GuideLiner catheter in use
- Y-adaptor with hemostasis valve (Tuohy-Borst type)
- Guidewire with diameter ≤ 0.014" / 0.36mm
- Sterile syringe (for flushing)
- Sterile heparinized saline (for flushing)

PREPARATIONS FOR USE

1. Prior to use, carefully inspect the GuideLiner packaging and components for damage.
2. Using sterile technique, transfer the dispenser coil with the GuideLiner catheter into the sterile field.
3. Remove the GuideLiner catheter from the dispenser coil and thoroughly flush the lumen of the GuideLiner catheter from the distal tip with sterile, heparinized saline solution.

DEPLOYMENT PROCEDURE

Deploy the GuideLiner catheter according to the following steps:

1. Secure the previously inserted guidewire and backload the distal tip of the GuideLiner catheter onto the guidewire and advance until the catheter is just proximal to the hemostasis valve.
2. Open the hemostasis valve and advance the GuideLiner catheter through the hemostasis valve and into the guide catheter.
3. Under fluoroscopy, advance the GuideLiner catheter beyond the distal tip of the guide catheter and into the desired location within the vessel.

WARNING: Never advance the GuideLiner catheter into a vessel with an effective diameter less than 2.5mm. Vessel injury, ischemia, and/or occlusion may result. If pressure in a vessel dampens after inserting the GuideLiner catheter, withdraw the catheter until the pressure returns to normal.

WARNING: Due to the size and non-tapered tip of the GuideLiner catheter, extreme care must be taken to avoid vessel occlusion and damage to the wall of the vessels through which this catheter passes.

4.

Using fluoroscopy, confirm the desired position of the GuideLiner catheter in the vessel.
5.

If performing an interventional procedure, backload the interventional device over the existing guidewire and advance the device through the guide catheter and GuideLiner catheter into the desired vascular space.
NOTE: If a second wire is used during the intervention and encounters resistance within the guide catheter, pull the wire back several centimeters and slowly readvance.
6.

Tighten the Y-adaptor hemostasis valve securely on the proximal shaft of the GuideLiner catheter to prevent back-bleeding.
7.

Perform the catheterization procedure. After completing the procedure, remove the GuideLiner catheter prior to removing the guide catheter from the vessel.
8.

Dispose of the GuideLiner catheter following standard hospital procedures.

STORAGE & HANDLING

No special storage or handling conditions.

LIMITED WARRANTY

Teleflex Medical LLC warrants that the GuideLiner catheter is free from defects in workmanship and materials prior to the stated expiration date. Liability under this warranty is limited to refund or replacement of any product, which has been found by Teleflex Medical LLC to be defective in workmanship or materials. Teleflex Medical LLC shall not be liable for any incidental, special or consequential damages arising from the use of the GuideLiner catheter. Damage to the product through misuse, alteration, improper storage or improper handling shall void this limited warranty.

No employee, agent or distributor of Teleflex Medical LLC has any authority to alter or amend this limited warranty in any respect. Any purported alteration or amendment shall not be enforceable against Teleflex Medical LLC.

THIS WARRANTY IS EXPRESSLY IN LIEU OF ALL OTHER WARRANTIES, EXPRESSED OR IMPLIED, INCLUDING ANY WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE OR ANY OTHER OBLIGATION OF TELEFLEX MEDICAL LLC.

PATENTS AND TRADEMARKS

May be covered by one or more U.S. or international patents.

See: www.teleflex.com/patents-intv

GuideLiner is a registered trademark of Teleflex Innovations S.à.r.l., or Teleflex Medical LLC, each a part of Teleflex Incorporated.

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