

SuperCross™ Microcatheters


Superb Support for Navigation of Tortuous Anatomy



Teleflex™
Empowering the future of healthcare

Superb Support in Challenging Cases

SuperCross Microcatheters provide support in tortuous anatomy and navigating bifurcated vessels.

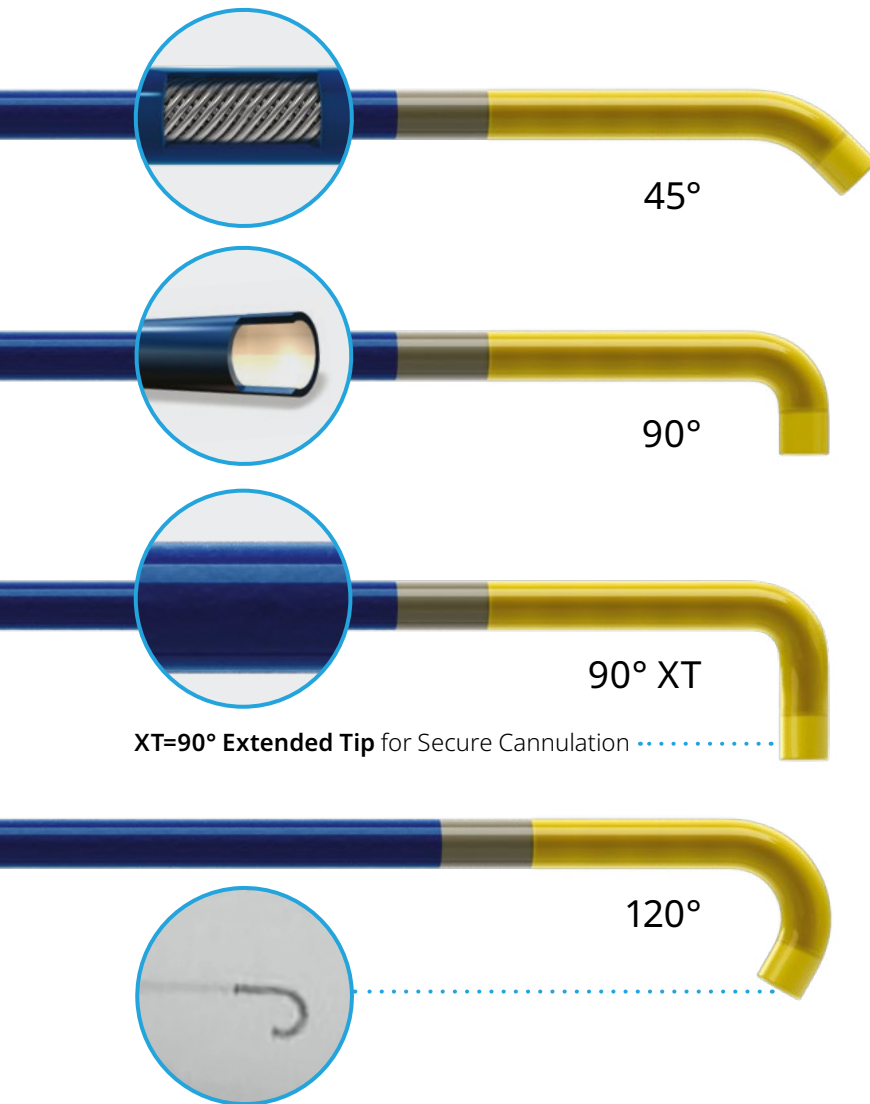

Angled Tip Catheters

Dual Coil Shaft Construction

Provides excellent torque response, flexibility, pushability and kink resistance

PTFE Inner Layer

Allows for outstanding guidewire movement and delivery

Hydrophilic Coating on the Distal 80 cm

Enables exceptional trackability around tight bends and tortuous anatomy

Embedded Platinum/Tungsten Coil

Provides enhanced visibility along the entire angled tip



Stainless Steel Braid for flexibility, pushability and kink resistance

Flexible Tip (FT) Catheters

Flexible Tip (FT) Version

For increased trackability through tortuous anatomy


ORDERING INFORMATION

SuperCross™ Microcatheters

The SuperCross™ Microcatheter is intended to be used in conjunction with steerable guidewires to access discrete regions of the coronary and/or peripheral vasculature. It may be used to facilitate placement and exchange of guidewires and other interventional devices and to subselectively infuse/deliver diagnostic and therapeutic agents.

MODEL	GUIDEWIRE COMPATIBILITY	DISTAL TIP	WORKING LENGTH	SHAFT CONSTRUCTION	DISTAL O.D.	PROXIMAL O.D.	DISTAL I.D.	PROXIMAL I.D.	HYDROPHILIC COATING
5302	0.014"	45°	130 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm
5303	0.014"	45°	150 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm
5304	0.014"	90°	130 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm
5305	0.014"	90°	150 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm
5306	0.014"	120°	130 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm
5307	0.014"	120°	150 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm
5308	0.014"	90° XT	130 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm
5309	0.014"	90° XT	150 cm	Coiled	2.4 Fr. 0.79 mm / 0.031"	3.2 Fr. 1.07 mm / 0.042"	0.017" (0.43 mm)	0.018" (0.46 mm)	Distal 80 cm

Packaged in quantities of 1 unit per box.

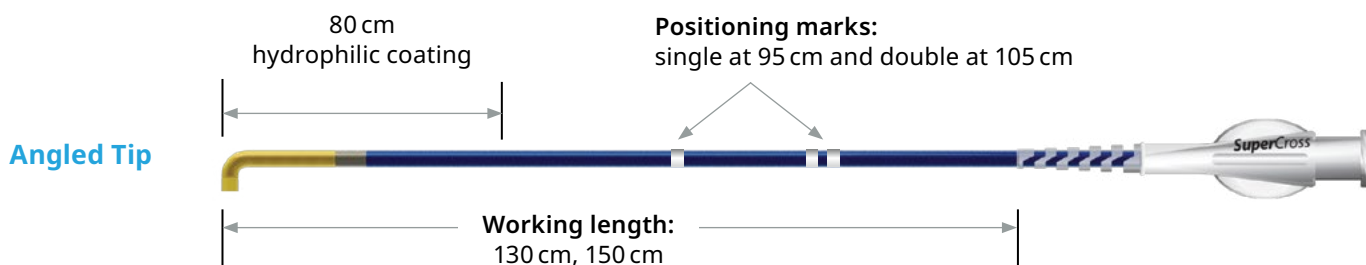


SuperCross™ Microcatheters—FT

The SuperCross™ Microcatheter is intended to be used in conjunction with steerable guidewires to access discrete regions of the coronary and/or peripheral vasculature. It may be used to facilitate placement and exchange of guidewires and other interventional devices and to subselectively infuse/deliver diagnostic and therapeutic agents.

MODEL	GUIDEWIRE COMPATIBILITY	DISTAL TIP	WORKING LENGTH	SHAFT CONSTRUCTION	DISTAL O.D.	PROXIMAL O.D.	DISTAL I.D.	PROXIMAL I.D.	HYDROPHILIC COATING
5340	0.014"	FT	130 cm	Braided	1.8 Fr. 0.61 mm / 0.024"	2.5 Fr. 0.84 mm / 0.033"	0.017" (0.43 mm)	0.021" (0.53 mm)	Distal 40 cm
5341	0.014"	FT	150 cm	Braided	1.8 Fr. 0.61 mm / 0.024"	2.5 Fr. 0.84 mm / 0.033"	0.017" (0.43 mm)	0.021" (0.53 mm)	Distal 40 cm

Packaged in quantities of 1 unit per box.

FT= Flexible Tip

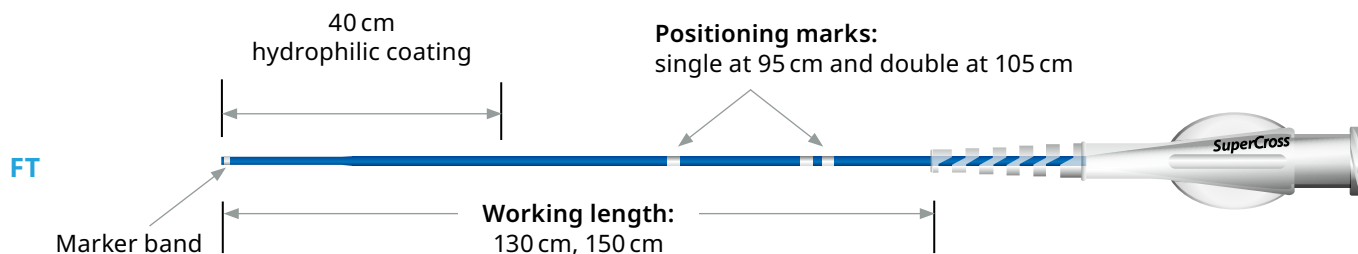

About Teleflex

As a global provider of medical technologies, Teleflex is driven by our purpose to improve the health and quality of people's lives. Through our vision to become the most trusted partner in healthcare, we offer a diverse portfolio with solutions in the therapy areas of anesthesia, emergency medicine, interventional cardiology and radiology, surgical, vascular access, and urology. We believe that the potential of great people, purpose driven innovation, and world-class products can shape the future direction of healthcare.

Teleflex is the home of Arrow™, Barrigel™, Deknatel™, LMA™, Pilling™, QuikClot™, Rüsch™, UroLift™ and Weck™ – trusted brands united by a common sense of purpose.

At Teleflex, we are empowering the future of healthcare. For more information, please visit teleflex.com.

Corporate Office

Phone +1 610 225 6800, 550 E. Swedesford Road, Suite 400, Wayne, PA 19087, USA

Regional Offices

United States: Phone +1 919 544 8000, Toll Free 866 246 6990, cs@teleflex.com, 3015 Carrington Mill Boulevard, Morrisville, NC 27560, USA

Latin America: Phone +1 919 433 4999, la.cs@teleflex.com, 3015 Carrington Mill Boulevard, Morrisville, NC 27560, USA

International: Phone +353 (0)9 06 46 08 00, orders.intl@teleflex.com, Teleflex Medical Europe Ltd., IDA Business and Technology Park, Dublin Road, Athlone, Co Westmeath, Ireland

CAUTION: Federal (USA) law restricts these devices to sale or use by or on the order of a physician.

Not all products are available in all regions. Please contact customer service to confirm availability in your region.

Teleflex, the Teleflex logo, and SuperCross are trademarks or registered trademarks of Teleflex Incorporated or its affiliates in the U.S. and/or other countries. Refer to the Instructions for Use for a complete listing of the indications, contraindications, warnings, and precautions. Information in this material is not a substitute for the product Instructions for Use. Not all products may be available in all countries. Please contact your local representative. For more information, please visit teleflex.com.

© 2025 Teleflex Incorporated. All rights reserved. MC-004044 Rev 0.1. Revised 10/2025.

Teleflex™
Empowering the future of healthcare